

Physicochemical Stability of diluted “Akynzeo®” Infusion Solutions in prefilled 0.9% Sodium Chloride Infusion Bags

¹Department of Pharmacy, University Medical Center of the Johannes Gutenberg University, Mainz, Germany

Zoe Mittmann¹, Judith Thiesen¹, Irene Krämer¹
zoe.mittmann@unimedizin-mainz.de

Background and Importance


Akynzeo® 235 mg/0.25 mg concentrate for solution for infusion, indicated for prevention of acute and delayed chemotherapy-induced nausea and vomiting, contains a combination of fosnetupitant (FOS) and palonosetron (PAL). Prior to administration, the Akynzeo® concentrate (20 mL vial) is diluted with either 30 mL or 100 mL 0.9% sodium chloride infusion solution. According to the SmPC, ready-to-administer (RTA) infusion solutions are physicochemically stable for 24 hours stored at room temperature [1].

Aim and Objectives

The aim of the study was to determine the **long term physicochemical stability** of two different ready-to-administer FOS/PAL infusion solutions in polyolefin (PO) infusion bags **stored at room temperature** over a **14-day period**.

Materials and Methods

Preparation of test solutions: FOS/PAL (Akynzeo® 235 mg/0.25 mg injection solution, 20 mL) diluted with 0.9% sodium chloride in prefilled PO bags

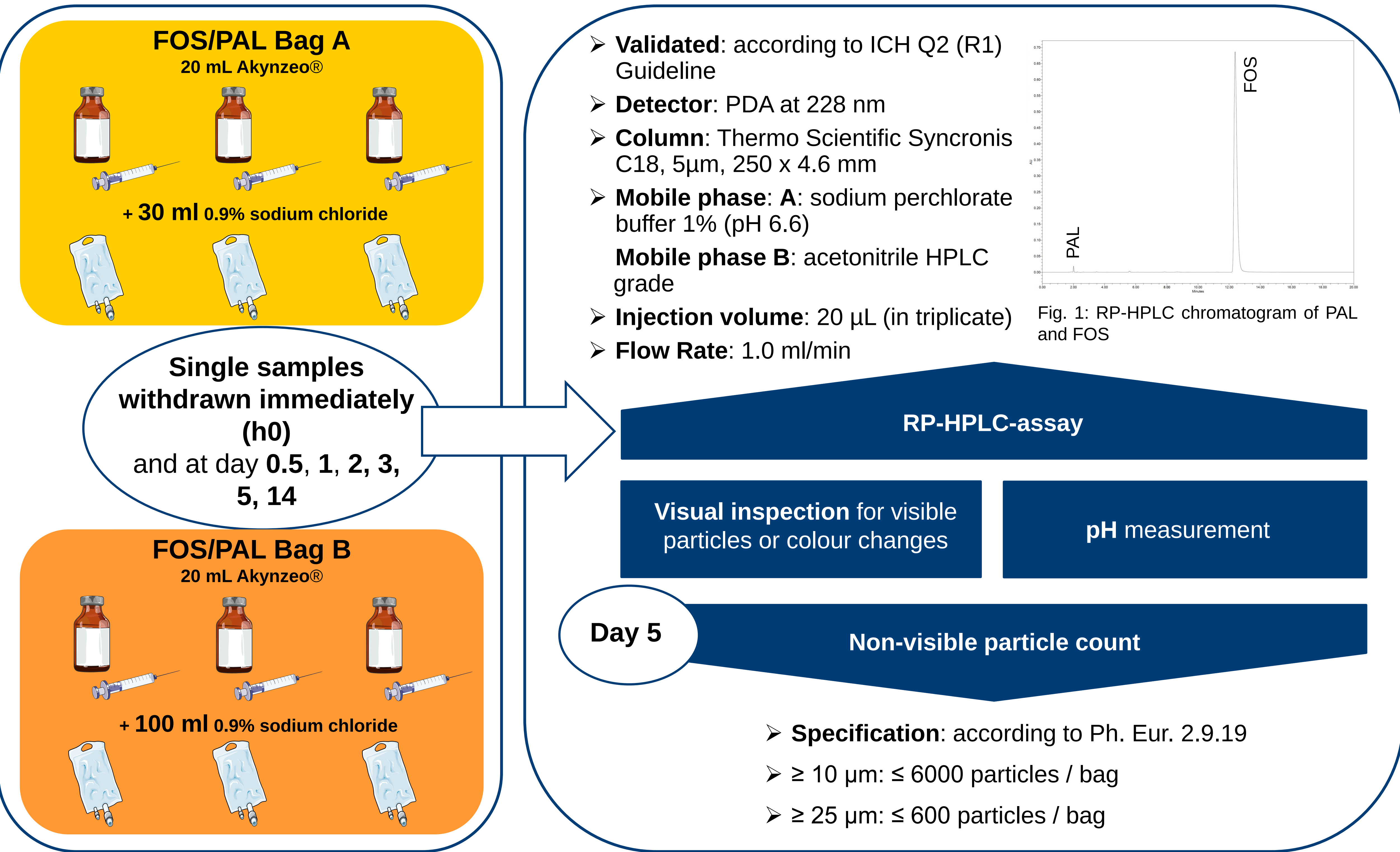

Results

Tab. 1: **FOS** concentration in **FOS/PAL test bag A**, stored at 25 °C over 14 days.

Nominal FOS concentration 4.70 mg/mL							
Test bag	Measured FOS concentration [mg/mL]						
	Day						
	0	0.5	1	2	3	5*	14
Bag A1	4.86	4.83	4.80	4.82	4.82	(4.62)	4.89
	4.88	4.83	4.80	4.83	4.83	(4.60)	4.91
	4.85	4.82	4.82	4.83	4.82	(4.60)	4.90
Bag A2	4.91	4.81	4.83	4.77	4.86	(4.62)	4.88
	4.90	4.81	4.83	4.77	4.86	(4.60)	4.88
	4.90	4.80	4.83	4.77	4.95	(4.59)	4.87
Bag A3	4.81	4.73	4.83	4.84	4.79	(4.47)	4.89
	4.81	4.73	4.83	4.83	4.79	(4.48)	4.88
	4.80	4.71	4.84	4.81	4.79	(4.45)	4.88
Mean [mg/mL] (n=9)	4.86	4.79	4.82	4.81	4.82	(4.56)	4.89
Percentage rate initial concentration [%]	100.00**	98.50	99.29	98.94	99.29	(93.83)	100.57
SD [mg/mL]	0.0404	0.0451	0.0130	0.0292	0.0281	(0.0659)	0.0116
RSD [%]	0.83	0.94	0.27	0.61	0.58	(1.45)	0.24

Tab. 2: **PAL** concentration in **FOS/PAL test bag A**, stored at 25 °C over 14 days.

Nominal PAL concentration 5.00 µg/mL							
Test bag	Measured PAL concentration [µg/mL]						
	Day						
	0	0.5	1	2	3	5	14
Bag A1	5.55	5.81	5.32	5.44	5.39	5.25	5.42
	5.55	5.78	5.34	5.46	5.45	5.37	5.42
	5.55	5.78	5.36	5.47	5.44	5.42	5.41
Bag A2	5.87	5.65	5.51	5.48	5.63	5.47	5.41
	5.85	5.67	5.53	5.51	5.63	5.48	5.44
	5.84	5.67	5.55	5.54	5.65	5.48	5.41
Bag A3	5.35	5.42	5.45	5.50	5.51	5.36	5.34
	5.34	5.43	5.46	5.51	5.50	5.38	5.32
	5.34	5.41	5.46	5.52	5.51	5.39	5.34
Mean [µg/mL] (n=9)	5.58	5.62	5.44	5.49	5.52	5.40	5.39
Percentage rate initial concentration [%]	100.00*	100.73	97.48	98.36	98.94	96.78	96.54
SD [µg/mL]	0.2090	0.1552	0.0798	0.0287	0.0879	0.0682	0.0419
RSD [%]	3.74	2.76	1.47	0.52	1.59	1.26	0.78

Tab. 3: pH values of **FOS/PAL test bag A**, stored at 25 °C over 14 days.

20 mL Akynzeo® + 30 mL 0.9% sodium chloride							
Test bag	pH value						
	Day						
	0	0.5	1	2	3	5	14
Bag A1	6.85	6.95	6.93	6.77	6.81	7.24	ND
Bag A2	7.10	7.01	7.11	7.01	6.97	7.17	ND
Bag A3	7.37	7.12	7.11	7.08	6.99	7.27	ND

Tab. 4: Non-visible particles of **FOS/PAL test bag A**, stored at 25 °C.

20 mL Akynzeo® + 30 mL 0.9% sodium chloride solution		
pH value		
Number [n] of non-visible particles on day 5		
Sample	≥ 10 µm	≥ 25 µm
Bag A1	6381*	32
Bag A2	4632	21
Bag A3	2240	16

Tab. 5: **FOS** concentration in **FOS/PAL test bag B**, stored at 25 °C over 14 days.

Nominal FOS concentration 1.96 mg/mL							
Test bag	Measured FOS concentration [mg/mL]						
	Day						
	0	0.5	1	2	3	5*	14
Bag B1	1.97	1.90	1.95	1.93	1.94	(1.79)	1.98
	1.98	1.90	1.95	1.93	1.94	(1.78)	1.98
	1.97	1.90	1.95	1.93	1.94	(1.77)	1.97
Bag B2	1.93	1.91	1.94	1.92	1.92	(1.77)	1.98
	1.93	1.91	1.94	1.92	1.91	(1.76)	1.97
	1.93	1.90	1.94	1.91	1.91	(1.74)	1.97
Bag B3	1.95	1.90	1.92	1.93	1.94	(1.74)	1.97
	1.95	1.89	1.93	1.93	1.94	(1.73)	1.97
	1.94	1.88	1.92	1.92	1.93	(1.72)	1.97
Mean [mg/mL] (n=9)	1.95	1.90	1.94	1.92	1.93	(1.75)	1.97
Percentage rate initial concentration [%]	100**	97.29	99.40	98.68	99.00	(89.97)	101.14
SD [mg/mL]	0.0179	0.0085	0.0115	0.0062	0.0106	(0.0216)	0.0035
RSD [%]	0.92	0.45	0.59	0.32	0.55	(1.23)	0.18

Tab. 6: **PAL** concentration in **FOS/PAL test bag B**, stored at 25 °C over 14 days.

Nominal PAL concentration 2.08 µg/mL							
Test bag	Measured PAL concentration [µg/mL]						
	Day						
	0	0.5	1	2	3	5	14
Bag B1	2.39	2.29	2.20	2.26	2.30	2.19	2.24
	2.36	2.29	2.20	2.25	2.30	2.19	2.24
	2.35	2.29	2.20	2.25	2.30	2.19	2.24
Bag B2	2.25	2.34	2.17	2.23	2.18	2.20	2.24
	2.21	2.33	2.18	2.22	2.18	2.20	2.24
	2.21	2.33	2.17	2.22	2.19	2.20	2.24
Bag B3	2.18	2.30	2.21	2.27	2.26	2.20	2.31
	2.18	2.28	2.21	2.27	2.25	2.21	2.31
	2.18	2.28	2.21	2.27	2.24	2.22	2.31
Mean [µg/mL] (n=9)	2.26	2.30	2.19	2.25	2.24	2.20	2.26
Percentage rate initial concentration [%]	100.00*	102.11	97.28	99.72	99.41	97.33	100.37
SD [µg/mL]	0.0805	0.0219	0.0164	0.0186	0.0472	0.0096	0.0342
RSD [%]	3.57	0.95	0.75	0.83	2.10	0.44	1.51

Tab. 7: pH values of **FOS/PAL test bag B**, stored at 25 °C over 14 days.

20 mL Akynzeo® + 100 mL 0.9% sodium chloride							
Test bag	pH value						
	Day						
	0	0.5	1	2	3	5	14
Bag B1	6.95	6.98	7.00	6.93	6.79	6.82	ND
Bag B2	6.85	6.96	6.97	6.90	6.73	6.78	ND
Bag B3	6.76	6.96	6.91	6.85	6.72	6.94	ND

Tab. 8: Non-visible particles of **FOS/PAL test bag B**, stored at 25 °C.

20 mL Akynzeo® + 100 mL 0.9% sodium chloride solution		
pH value		
Number [n] of non-visible particles on day 5		
Sample	≥ 10 µm	≥ 25 µm
Bag B1	3820	67
Bag B2	2253	13
Bag B3	2753	7

➤ **Visual inspection:** solutions of FOS/PAL test bags A and test bags B remained clear, without change of color or appearance of visible particles over 14 days.

Conclusion

Ready-to-administer FOS/PAL infusion solutions (Akynzeo® 235 mg/0.25 mg concentrate for solution for infusion) diluted with 30 mL and 100 mL in prefilled 0.9 % sodium chloride PO bags are physicochemically stable up to **14 days when stored at 25 °C** and protected from light.

References

[1] Akynzeo Summary of product characteristics. https://www.ema.europa.eu/en/documents/product-information/akynzeo-epar-product-information_en.pdf